



smart optics  
Lise-Meitner-Allee 10  
44801 Bochum  
Germany

☎ 0049 234 298280  
☎ 0049 234 2982820  
💻 [www.smartoptics.de](http://www.smartoptics.de)  
✉ [info@smartoptics.de](mailto:info@smartoptics.de)

# smart optics

## Activity 885 Mark 2

### Fast

In just 29 seconds, the  
Activity 885 Mark 2  
scans a complete jaw

### Efficient

Articulated  
scanning with  
Activity 885 Mark 2

### Uncomplicated

Output in an open STL  
and PLY file format  
and simple handling



✓ dental Scan

✓ secondDie

✓ Texture scan monochrome

🔍 multiDie

🔍 multiCase

🔍 Fixator scan

🔍 Triple Tray® impression scan

Compatible with exocad® Dental CAD

# Activity 885 Mark 2


## Technical Data

Dimension (B x H x T)	490 x 430 x 440 mm
Measuring field (X x Y x Z)	80 x 60 x 85 mm
Resolution	1,3 MP
Accuracy (according to ISO 12836)	6 µm
Weight	32 kg
Power Supply Voltage	100 - 240 VAC, 50/60 Hz
Interfaces	2x USB
Housing Color	anthracite
Scope of delivery	1x Activity 885 Mark 2 1x Object holder 1x 3D calibrations model 2x USB cable 1x Mains cable 1x Data carrier

## Features

dental Scan	✓
secondDie	✓
Texture scan <small>monochrome</small>	✓
Compatible with exocad® Dental CAD	✓
Fixator scan	🔑
Triple Tray® impression scan	🔑
multiCase	🔑
multiDie	🔑

Be smart. Be professional.

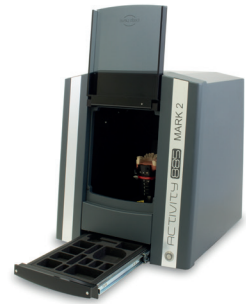

## Measurements

### Complete jaw

Scanning	16 sec
Postprocessing	13 sec
Total	29 sec

### Single tooth

Scanning	27 sec
Postprocessing	9 sec
Total	36 sec

### 3-unit bridge

Scanning	38 sec
Postprocessing	22 sec
Total	60 sec

## Contact us for further information



smart optics  
Lise-Meitner-Allee 10  
44801 Bochum  
Germany

☎ 0049 234 298280  
☎ 0049 234 2982820  
🌐 [www.smartoptics.de](http://www.smartoptics.de)  
✉ [info@smartoptics.de](mailto:info@smartoptics.de)

📘 smartopticsDE  
🐦 smartopticsDE  
📺 SmartOpticsVideos

# Activity 885 Mark 2

Be smart. Be professional.

For all-rounders:  
Activity 885 Mark 2



Not every successor is successful, well, this one is. The fully automated Activity 885 Mark 2 is the successor model of the popular Activity 885. The two models differ significantly, especially in terms of speed. The Activity 885 Mark 2 is up to 50 percent faster than its predecessor. Giving you extra time to allow increased workload and more efficiency.

The scanner also scores with its comfortable fixtures which also include the secondDie software module next to a 36 month guarantee. With secondDie, the scanning process is a further 40 percent faster. The Activity 885 Mark 2 completes the data for a segmented 3-unit bridge in only 40 seconds.

Of course the Activity 885 Mark 2 can scan articulated jaw models in a condyle-related manner using the smart optics scan fixators and the data to the dental CAD.

The Activity 885 Mark 2 captivates with its high quality, both in terms of workmanship as well as with its scanning results.



- **Fast**  
In just 29 seconds, the Activity 885 Mark 2 scans a complete jaw
- **Uncomplicated**  
Output in an open STL and PLY file format and simple handling
- **Efficient**  
Articulated scanning with Activity 885 Mark 2

Contact us for further information



smart optics  
Lise-Meitner-Allee 10  
44801 Bochum  
Germany

☎ 0049 234 298280  
☎ 0049 234 2982820  
💻 [www.smartoptics.de](http://www.smartoptics.de)  
✉ [info@smartoptics.de](mailto:info@smartoptics.de)

📘 smartopticsDE  
📺 smartopticsDE  
📺 SmartOpticsVideos